

mectron

→ MECTRON EASYJET PRO







# EASYJET PRO

## airpolishing for turbine coupling

→ easyjet pro opens up a new dimension of turbine couplings. Just a simple click and your turbine becomes a skillful prophylaxis unit, which offers the same well-proven features as mectron's table top units: a powerful powder jet, easy handling and reliable airpolishing. Shape and balance also enable a good handling and a controlled precise powder jet direction.

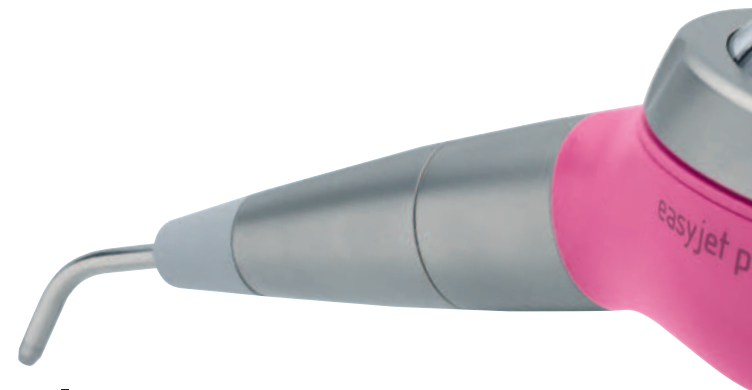

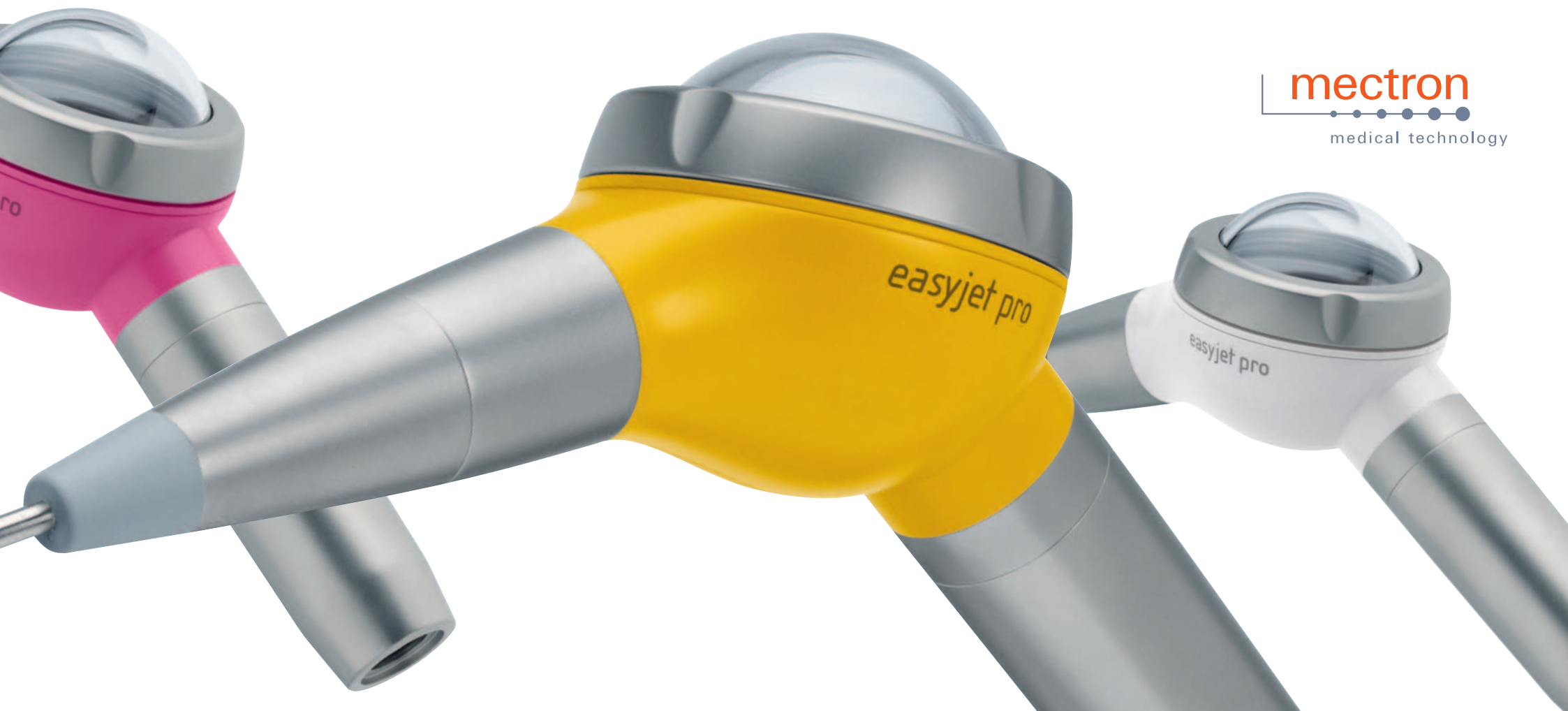

→ COLORS  
The easyjet pro is available in seven different colors. The colorful selection will become a prophylaxis highlight in your office!

# → THE TECHNOLOGY

reliable and ergonomic

→ easyjet pro stands out due to its innovative technology, which allows high reliability. All parts are optimized for daily use.



## → SPRAY HANDPIECE

- 120° angle
- just a click to connect
- 360° rotating head
- powder and water travel separately until they reach the nozzle
- the water surrounds the jet of powder and air and prevents powder dust
- fully sterilizable



## → POWDER CHAMBER

- powder capacity of 14 grams
- a new technique allows a constantly effective powder jet
- no clogging

## → TECHNICAL DATA

- incoming air pressure between 2.3–3.2 bar
- incoming water pressure between 0–3 bar

CE  
0476



- SILICON PROTECTOR
- protects teeth
  - prevents any pollutant entry
  - dismantable
  - fully sterilizable



- POWDER CHAMBER CAP
- stainless steel ring-nut
  - transparent dome for optimal control of powder level
  - triangular thread for easy opening
  - thanks to its specific shape there will be no powder deposits on the thread





# THE TURBINE ADAPTER

→ easyjet pro offers you a high degree of flexibility with its six different turbine coupling connectors. With just a click you can carry out the prophylaxis treatment. The connector can easily be changed by an authorized mectron service in case of a turbine coupling change.



→ BIEN AIR UNIFIX  
AND UNIFIX L



→ KAVO MULTIFLEX  
AND MULTIFLEX LUX



→ SIRONA R AND B



→ W & H ROTO QUICK  
AND ROTO QUICK LUX



→ MIDWEST



→ BORDEN

# → THE POWDER

→ mectron prophylaxis powder is delivered in convenient plastic bottles. Each bottle contains 250 grams of powder. The package includes two bottles of mectron prophylaxis powder.



## → EFFECT

Particles of sodium bicarbonate at a high rate of acceleration combined with warm water enable an extremely fine, regular and precise jet. The jet eliminates discolorations more effectively, gently and faster than when using rubber heads and polishing pasta, even in the interdental spaces.



Before, during and after easyjet pro-treatment

## → QUALITY

Its optimized flowing and hydrophobic feature avoid clogging and allow a more regular powder jet with amazing treatment results.



## → TASTE

The new taste is originated by a new lemon flavor mixed with a larger quantity of sweet saccharin.



## → POWDER DOSAGE

Each bottle has its dosage measuring cap, indicating the proper powder quantity for each mectron airpolisher.

# MECTRON EASYJET PRO PLUG AND SPRAY!

[www.mectron.com](http://www.mectron.com) or [mectron@mectron.com](mailto:mectron@mectron.com)

mectron s.p.a., via Loreto 15/A, 16042 Carasco (Ge), ITALIA, tel +39 0185 35361, fax +39 0185 351374



**mectron**  
medical technology

© Copyright mectron S.p.A., Carasco, Italy  
All rights reserved. Texts, pictures and graphics of mectron brochures are protected by copyright and other protection laws. Without written approval of mectron S.p.A. the contents may not be copied, distributed, changed or made available to third parties for commercial purposes.

DEP1013EN1105